

# DENERVEX™

ROTATIONAL ABLATION DENERVATION DEVICE

## ENDURING RELIEF™



# Introducing the DenerveX™ System

A procedural approach designed for the long-lasting relief of facet joint pain.

Approximately 31% of the chronic low back pain patients report is attributed to the facet joint.<sup>1</sup> DenerveX™ is a new treatment option for facet joint pain.


The DenerveX™ treatment uses Rotacapsulation™, a combination of high heat and rotational capsular tissue shaving, in a minimally invasive posterior procedure for the treatment of Facet Joint Syndrome.

## + LONG-LASTING PAIN RELIEF AND IMPROVED MOBILITY<sup>2</sup>

The DenerveX™ design is based on a manual procedure with precision and reproducibility built-in.<sup>2</sup>


- ▶ The patient pain score, as measured by the VAS pain scale, at 3 years follow up was reduced 75% to 100% in more than 50% of patients treated for cervical, thoracic, and/or lumbar facet joint pain.
- ▶ The patient evaluated disability score, as measured by the Oswestry Disability Index, at 3 years follow up was improved 75% to 100% in more than 50% of patients treated for facet joint pain.

Percent of patients experiencing a  $\geq 50\%$  improvement in pain at 3-year follow up<sup>1</sup>

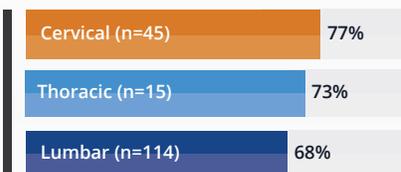

The majority of patients achieved a 75% to 100% improvement (VAS pain score) at a minimum 3-year follow up.

Percent of patients experiencing a  $\geq 50\%$  improvement in mobility at 3-year follow up<sup>1</sup>

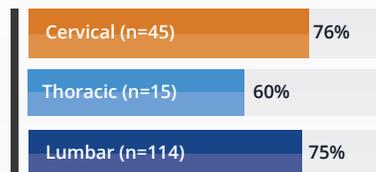

Mobility also improved. Most patients experienced a 75% to 100% improvement in the Oswestry Disability Index at long term follow up.

## + EASY-TO-USE


The physician designed and tested DenerveX™ system incorporates procedural simplicity, including:

- ▶ Minimally invasive procedure
- ▶ Three Effects, based on patients BMI, designed to deliver precise therapy for each patient
- ▶ Automatic control of therapy time
- ▶ All inclusive, ready-to-use, single use device kit

DenerveX™ Pro-40 Generator



DenerveX™ Procedure Kit

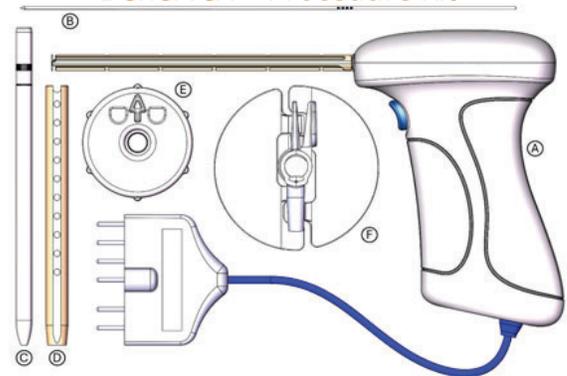

- ▶ Hand Piece (A), K-Wire (B), Procedure Dilator (C)
- ▶ Portal Tube (D), Portal Driver (E), Tissue Stabilizer (F)

## + INNOVATIVE TREATMENT FOR FACET JOINT SYNDROME

The DenerveX™ system is exclusively designed as a system to treat facet joints. It uses Rotacapsulation™ to deliver high heat and capsular tissue shaving in a posterior capsulectomy procedure that:

- ▶ Treats the posterior synovial capsule
- ▶ Removes the facet joint nerve receptors
- ▶ Terminates the nerve





## About Us

Medovex was formed to acquire and develop a diversified portfolio of potentially ground breaking medical technology products. Criteria for selection include those products with potential for significant improvement in the quality of patient care combined with cost effectiveness. The Company's first pipeline product, the DenerveX™ system, is intended to provide long-lasting relief from pain associated with facet joint syndrome at significantly less cost than currently available options.

## Contact Us

**Medovex**

**Unit 1/18 Dequetteville Terrace**

**Kent Town**

**SA 50670**

**Adelaide**

**Australia**

**[medovex.com](http://medovex.com)**

### References:

1. *Laplante B L, Ketchum J M, Saullo T R and DePalma M J. Multivariable Analysis of the Relationship Between Pain Referral Patterns and the Source of Chronic Low Back Pain. Pain Physician 2012, 15: 171-178.*
2. *Haufe, S M, Mork A R. Endoscopic Facet Debridement for the treatment of facet arthritic pain – a novel new technique. Int J Med Sci 2010; 7(3): 120 – 123.*



DenerveX™ is currently not CE Marked or commercially available in the United States of America.

*Medovex and DenerveX are trademarks of Medovex Pty Ltd.*

© Copyright 2025 Medovex Pty Ltd All rights reserved.